Autoclave Monitor Tape

.75 in x 60 yd / 19 mm x 54 m

EN

INSTRUCTIONS FOR USE

Product Description

The HuFriedyGroup Autoclave Monitor Tape is designed for use with steam sterilization as a method of securing packs and identifying processed from unprocessed items.

Intended Use

The HuFriedyGroup Autoclave Monitor Tape is designed to demonstrate that the pack has been exposed to a sterilization process and to distinguish between processed and unprocessed packs. After steam sterilization, the chemical indicator lines will show a visual color change from off yellow to dark brown/black.

Use the HuFriedyGroup Autoclave Monitor Tape with the following steam sterilization cycles.

- 1. 121°C, 30 Minutes; Gravity Steam
- 2. 132°C, 4 Minutes; Dynamic Air Removal Steam
- 3. 135°C, 3 Minutes; Dynamic Air Removal Steam

Instructions for Use

When using the HuFriedyGroup Autoclave Monitor Tape to secure sterilization wraps it is recommended to:

- 1. Remove the product label and plastic packaging and install the tape roll in a compatible tape dispenser (HuFriedyGroup Tape Dispensers IMS-1270 or IMS-1270H).
- 2. Unwind a section of Autoclave Monitor Tape long enough to secure the pack.
- 3. Place one end of the Autoclave Monitor Tape on the package and apply finger pressure across the length of the tape to secure the pack. Do not stretch the tape.
- 4. Repeat steps 2 and 3 if the method of pack wrapping / size requires two or more strips of indicator tape for secure closure. Do not allow any openings in the folds.
- 5. When the steam sterilization cycle is complete, the chemical indicator lines will show a visual color change from yellow to dark brown/black.
- 6. When the sterile pack is ready for use, open by using aseptic technique to avoid contaminating the internal contents.


Shelf-Life:

Unused: 24 Months from date of manufacture.

Available Options:

Part Code	Procedure Color
IMS-1253	Cement
IMS-1254	Composite
IMS-1255	Endo
IMS-1256	Exam
IMS-1257	Extraction
IMS-1258	Perio
IMS-1259	Prophy
IMS-1260	Pedo
IMS-1262	Gray
IMS-1263	Orange
IMS-1264	Purple

Part Code	Procedure Color
IMS-1265	Yellow
IMS-1266	Red
IMS-1267	Green
IMS-1268	Blue
IMS-1269	Black
IMS-1243	Orthodontic
IMS-1247	Restorative
IMS-1248	Surgery
IMS-1251	Amalgam
IMS-1252	C and B

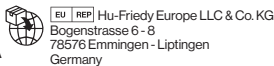Storage Conditions

59°F/5°C — 86°F/30°C

ISO 11140-1 Type 1

<50%

Do not freeze.
Keep away from sterilants and excessive heat.

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© 2026 Hu-Friedy Mfg. Co., LLC. Tous droits réservés.
© 2026 Hu-Friedy Mfg. Co., LLC. Alle rechten zijn voorbehouden.