Intended use: Designed to verify the cleaning performance in the tank of ultrasonic irrigators and cleaners.

Instructions for Use:


Test to be conducted daily, in an EMPTY tank.

1. Ensure that the Ultrasonic Test Holders **REF** U505 are clean and dry.
2. Ensure gloves are clean and dry before handling the indicator.
3. Place 1 VERIFY® Ultrasonic Indicator in each of the 5 recommended holders.
The indicator should be placed in the middle of the spring inside of the holder.


4. We recommend attaching 1 holder to each of the 4 corners of the tray using the clip on the back of the holder and recommend setting 1 holder in the center of the tray. Example of recommended arrangement in tray shown below:

Medium to Large Ultrasonic Trays

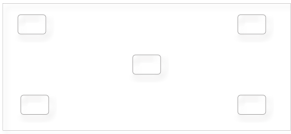Small Ultrasonic Trays

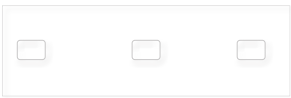5. Place tray into the ultrasonic cleaner or sonic irrigator and run a cycle with at least 5 minutes of sonic activity. For sonics with multiple levels, place a total of five holders as shown, across all levels.
6. After the cycle is completed, remove the holders and use clean gloves to remove the Indicator from each holder. Caution should be used, as the device may be hot and any residual soil may stain.

7. Inspect the indicator for evidence of soil by placing it against a white background.

Single use only, do not reuse processed/partially processed indicators due to a risk of a false Pass result.

DISPOSAL: This product may be disposed of as regular waste.

Refer to the STERIS symbol glossary for definitions at <https://www.steris.com/symbol-glossary>.

Note: Serious incidents that have occurred in relation to this medical device should be reported to the manufacturer and competent authority in the country where the incident occurred.

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