


VERIFY™ ULTRASONIC INDICATOR

Ultrasonic Monitoring



NO CUTTING CORNERS WITH YOUR ULTRASONIC CYCLE

Advancements in surgical instruments have dramatically improved patient outcomes. With increased reliance on ultrasonic cleaners to keep pace with caseloads, ultrasonic monitoring has never been more important.

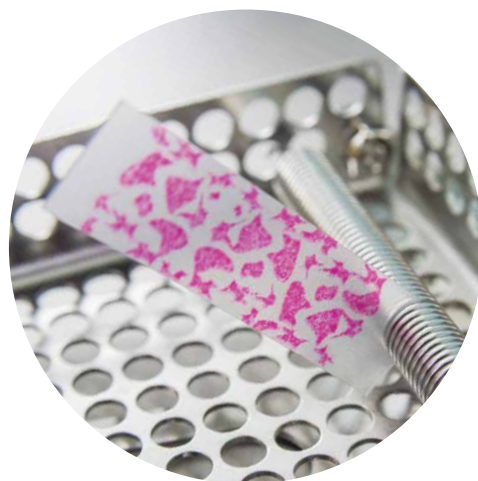

COMPLETE ASSURANCE

Ultrasonic cleaning uses cavitation, chemistry, time and temperature to remove soils after manual cleaning. Problems such as insufficient sonic energy, water level, improper temperature and degassing may have an impact on cycle outcomes.

THE VERIFY ULTRASONIC INDICATOR DETECTS THE FOUR KEY PARAMETERS OF ULTRASONIC CLEANING:

Zone Monitoring

helps you measure Ultrasonic Activity throughout the tank and not just one area

Detects Chemistry Presence

to ensure active soil removal

Monitors Time

for successful soil removal and monitors temperature for ideal enzymatic performance

Drives Compliance

and quality assurance by monitoring the performance of your Ultrasonic Cleaner or Sonic Irrigator



Realistic Challenge

is produced by the holder's design that mimics intricate surgical instruments



Safe to Handle

synthetic test soil mimics blood and common tissues found on devices

ORDERING INFORMATION

PART NUMBER	DESCRIPTION	QUANTITY
U504	VERIFY Ultrasonic Indicators	200/pack
U505	VERIFY Ultrasonic Indicator Holder	5/pack

ANSI/AAMI ST79: 13.2, 2017
Mechanical cleaning equipment should be tested upon installation, each day that it is used, and after major repairs.

TO LEARN MORE ABOUT THE VERIFY™ ULTRASONIC INDICATOR, please contact your local STERIS representative.

healthcare.steris.com
Mentor, OH 44060 | USA

